



# First MPJ Hemi Implant

The BioPro First MPJ Hemi Implant is the globally trusted **pre-fusion option** for the treatment of hallux rigidus/limitus. BioPro offers the most comprehensive phalangeal based hemi implant system, along with the addition of our new patented Hemi**EDGE**.™

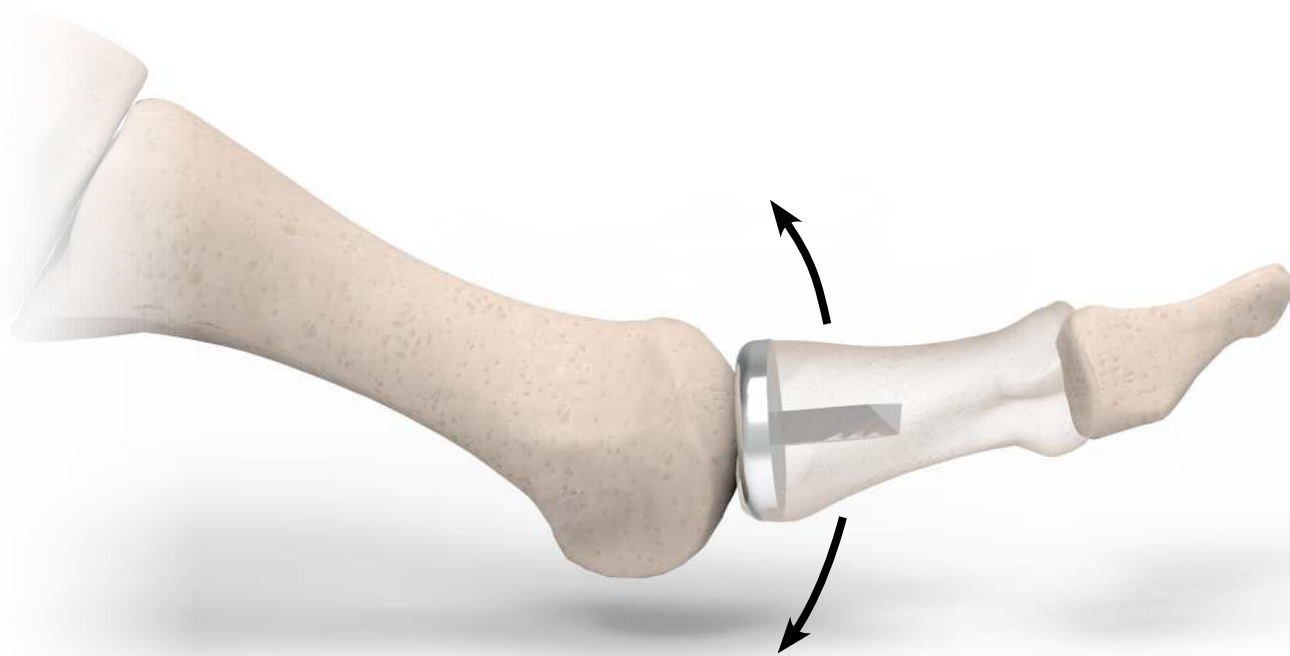

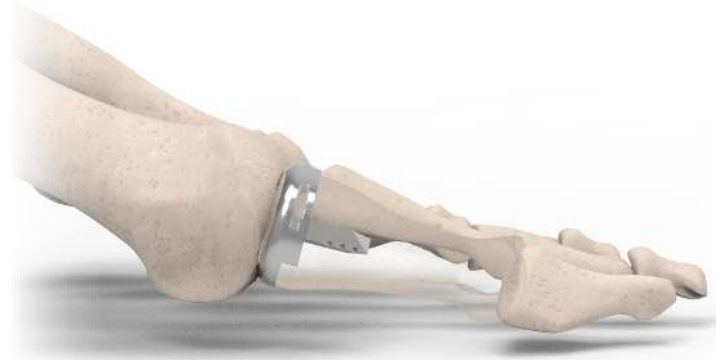

## Reproducible results

- On average studies report 95%+ implant survivorship\*
- 65+ years of successful, clinical use<sup>9</sup>
- Simple instrumentation

## Satisfied patients

- 97% patient satisfaction at 8.3 years<sup>1</sup>
- Immediate weight bearing at the discretion of the surgeon
- Return to activity on average 5 weeks faster than fusion patients<sup>1</sup>

Potential complications of hemiarthroplasty include implant subsidence, bone overgrowth and implant migration. The HemiEDGE was designed to address these complications by encompassing the proximal phalanx with an overlapping edge around the medial, lateral and dorsal sides. The implant articulating surface remains consistent with the standard design.

## A trusted pre-fusion option

In the right patient, our implants have been proven to last **20+ years**<sup>9</sup>. If complications do arise, studies show conversion from a failed hallux phalangeal hemiarthroplasty to arthrodesis have similar success rates to primary arthrodesis.<sup>13</sup>

System Feature	Benefit
Low Profile	Minimum bone removal doesn't sacrifice FHB
Barbed Stem	Secure fixation that doesn't require drilling or bone cement
Multiple Variations	Porous Coated (PC) for improved biological fixation and Titanium for metal sensitivities

	17MM	18.5MM	20MM	21.5MM	23MM
<b>Cobalt Chrome NPC</b>	10060	17033	10061	14958	10062
<b>Cobalt Chrome PC</b>	10412	17034	10413	14960	10414
<b>Titanium NPC</b>	16813	17198	16814	16815	16816
<b>Titanium PC</b>	17035	17197	17036	17037	17038
<b>HemiEDGE</b>	19538	19539	19540	19541	19542

Visit [www.bioproimplants.com](http://www.bioproimplants.com) or call us at (810) 982-7777

1. Timothy Voskuijl, MD, Ron Onstenk, MD. Operative Treatment for Osteoarthritis of the First Metatarsophalangeal Joint: Arthrodesis Versus Hemiarthroplasty. The Journal of Foot & Ankle Surgery xxx (2015) 1-4
2. Karin H. Simons, MD, Pieter van der Woude, MD, Frank W.M. Faber, MD, PhD, Paulien M. van Kampen, PhD, Bregje J.W. Thomassen, PhD. Short-Term Clinical Outcome of Hemiarthroplasty Versus Arthrodesis for End-Stage Hallux Rigidus. The Journal of Foot & Ankle Surgery xxx (2015) 1-4
3. Giza E, Sullivan MR. First Metatarsophalangeal Hemiarthroplasty for Grade III and IV Hallux Rigidus Techniques in Foot and Ankle Surgery 4(1):10-17,2005
4. Clement, N. D., MacDonald, D., Dall, G. F., Ahmed, I., Duckworth, A. D., Shalaby, H. S., & McKinley, J. (2016). Metallic hemiarthroplasty for the treatment of end-stage hallux rigidus. Bone Joint J, 98-B(7), 945-951.
5. Taranow, DO, et al. Contemporary Approaches to Stage II and III Hallux Rigidus: The Role of Metallic Hemiarthroplasty of the Proximal Phalanx. Foot and Ankle Clinics, Volume 10, Issue 4, 713 - 728
6. Roukis TS, Townley, MD. BIOPRO resurfacing endoprosthesis versus periarticular osteotomy for hallux rigidus: short-term follow-up and analysis. Journal of Foot & Ankle Surgery 2003;42(6):350-8
7. Juan C. Goetz, DPM, Charles O. Townley MD, Warren Taranow, DO. An Update on the Metallic Hemiarthroplasty Resurfacing Prosthesis for the Hallux. Presented at the 56th Annual Meeting and Scientific Seminar of the American College of Foot and Ankle Surgeons. Orlando FL February 1998
8. Charles G. Kissel, DPM, FACFAS, Zeeshan S. Husain, DPM AACFAS, Paul H. Wooley, PhD, Michael Kruger, MS, Mark A. Schumaker, DPM, Michael Sullivan, DPM, and Todd Snoeyink, DPM. A Prospective Investigation of the Biopro® Hemi-Arthroplasty for the First Metatarsophalangeal Joint. The Journal of Foot & Ankle Surgery 47(6):505-509, 2008
9. Taranow, DO, Townley, MD. Metallic proximal phalangeal hemiarthroplasty for hallux rigidus. Operative Techniques in Orthopaedics 1999;9(1):33
10. Townley, MD, Taranow, DO. A metallic hemiarthroplasty resurfacing prosthesis for the hallux metatarsophalangeal joint. Foot & Ankle International 1994;15(11):575-80
11. Christine C. Salonga, DPM, David C. Novicki, DPM, FACFAS, Martin M. Pressman, DPM, FACFAS, D. Scot Malay, DPM, MSCE, FACFAS. A Retrospective Cohort Study of the BioPro Hemiarthroplasty Prosthesis. The Journal of Foot & Ankle Surgery 49 (2010) 331-339
12. Giza, E., Sullivan, M., Ocel, D., Lundeen, G., Mitchell, M., & Frizzell, L. (2010). First metatarsophalangeal hemiarthroplasty for hallux rigidus. International Orthopaedics, 34(8), 1193-1198. <http://doi.org/10.1007/s00264-010-1012-x>
13. Garras, D. N., Durinka, J. B., Berck, M., Miller, A. G., & Raikin, S. M. (2013). Conversion Arthrodesis for Failed First Metatarsophalangeal Joint Hemiarthroplasty. Foot & Ankle International, 34(9), 1227-1232.

\* Implant survivorship average from studies 1-12

HemiEDGE Patent No. 8,845,750



+1-810-982-7777  
info@bioproimplants.com  
2929 Lapeer Road  
Port Huron, MI 48060



CE Standard Implants Only

MKT15 01